


**DISPOSABLE INTRAVENOUS
CATHETER**

**INFUSION / TRANSFUSION
THERAPY**

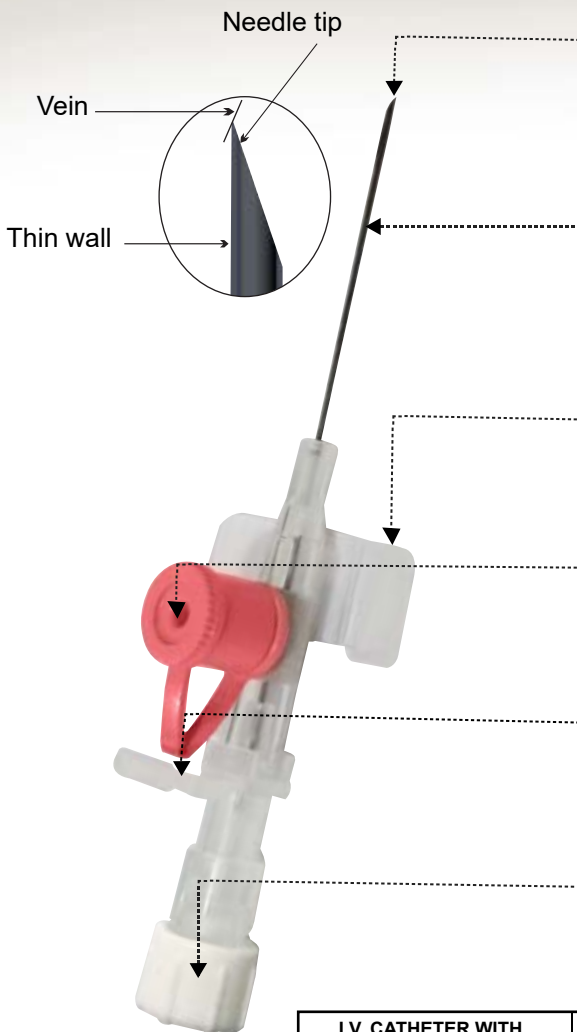

Sharp 3-Facet Back-Cut Needle

- Sharp and thin walled needle with reduced bevel angle for less contact area during insertion resulting **atraumatic cannulation**.

Fluorinated Ethylene Propylene (FEP) Catheter

- Special tapered and kink resistant **Fluorinated Ethylene Propylene (FEP)** radiopaque catheter for easy insertion with optimal flow rate.

Flexible Wing

- Special flexible wing to **provide easy and secure fixation**.

Injection Port

- Injection port with non-return Silicone Valve for **intermittent medication**.

Needle Grip

- Ergonomic needle grip design for **one-handed puncture technique** and comfortable handling.

Removable Closing Cone

- Luer lock plug at the end of catheter to **prevent blood contamination**.

- ✓ Disposable.
- ✓ Sterile.
- ✓ Easy to use.
- ✓ Latex free.
- ✓ PVC free.
- ✓ Individually packed in hard blister packaging.
- ✓ Packaging: 50pcs/box, 10 boxes / carton.

I.V. CATHETER WITH INJECTION PORT & WINGS	CATALOGUE NO.	SIZE	COLOR CODE	OD (mm)	ID (mm)	LENGTH (mm)	FLOW RATE (ml/min)
	IC2001-14	14G	●	2.0	1.7	50	305
	IC2001-16	16G	●	1.7	1.3	50	200
	IC2001-18	18G	●	1.3	0.9	45	95
	IC2001-20	20G	●	1.1	0.8	33	65
	IC2001-22	22G	●	0.9	0.6	25	36
	IC2001-24	24G	●	0.7	0.5	19	23
	IC2001-26	26G	●	0.6	0.45	19	17



Ideal Healthcare Sdn. Bhd.(484861-V)
No. 70 & 71, Jalan Sungai Tukang 2/1, Kawasan Perusahaan Sungai Tukang,
08000 Sungai Petani, Kedah, Malaysia.
Tel : 6 (04) - 422 7775 / 422 8775 Fax : 6 (04) - 422 6775
Email : info@idealcare.com.my

E. L. NO.: MDA-0003-KD314

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= ISO 13485 =

www.idealcare.com.my




idealcareMY

IHSB-CAT-BDD-ITT-13

**DISPOSABLE INTRAVENOUS
CATHETER**

INFUSION / TRANSFUSION
THERAPY

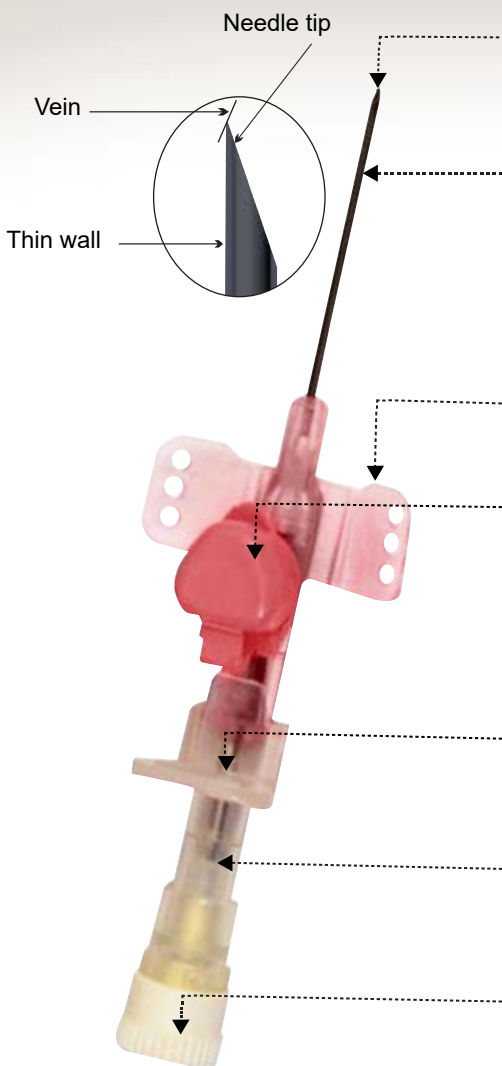

Sharp 3-Facet Back-Cut Needle

- Sharp and thin walled needle with reduced bevel angle for less contact area during insertion resulting **atraumatic cannulation**.

Biocompatible Polyurethane (PUR) Catheter

- Special tapered bio-compatible **Polyurethane (PUR)** catheter reduces the chances of phlebitis and thrombo-phlebitis; kink resistant allow longer indwelling time.
- Smooth surface catheter embedded with **radiopaque lines** for better x-ray visualization.

Flexible Perforated Colored Wing For Ventilation

- Special flexible perforated colored wing for **secure fixation and ventilation**.

User Friendly Snap Flip Cap

- Injection port with snap flip cap for **one-handed operation of easy opening and closing**.
- Compatible with both luer syringe hub.
- Needle-free injection reduces risk of needlestick injury.
- Universal color-coded on port for easy gauge size identification.

Needle Grip

- Ergonomic needle grip design for **one-handed puncture technique** and comfortable handling.

Transparent Flashback Chamber and Hydrophobic Filter Membrane

- Transparent flashback chamber for **easy visualization** of blood.
- Hydrophobic filter to prevent blood leakage during cannulation.

Removable Closing Cone

- Luer lock plug at the end of catheter to **prevent blood contamination**.

- ✓ Disposable.
- ✓ Sterile.
- ✓ Easy to use.
- ✓ Latex free.
- ✓ PVC free.
- ✓ Individually packed in hard blister packaging.
- ✓ Packaging: 50pcs/box, 10 boxes / carton.

I.V. CATHETER WITH INJECTION VALVE	CATALOGUE NO.	SIZE	COLOR CODE	OD (mm)	ID (mm)	LENGTH (mm)	FLOW RATE (ml/min)
	IC2002-14	14G	Orange	2.0	1.7	50	305
	IC2002-16	16G	Grey	1.7	1.3	50	200
	IC2002-18	18G	Green	1.3	0.9	45	95
	IC2002-20	20G	Red	1.1	0.8	33	65
	IC2002-22	22G	Blue	0.9	0.6	25	36
	IC2002-24	24G	Yellow	0.7	0.5	19	23
	IC2002-26	26G	Purple	0.6	0.45	19	17



Ideal Healthcare Sdn. Bhd.(484861-V)
No. 70 & 71, Jalan Sungai Tukang 2/1, Kawasan Perusahaan Sungai Tukang,
08000 Sungai Petani, Kedah, Malaysia.
Tel : 6 (04) - 422 7775 / 422 8775 Fax : 6 (04) - 422 6775
Email : info@idealcare.com.my

E. L. NO.: MDA-0003-KD314

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= ISO 13485 =

www.idealcare.com.my




idealcareMY

IHSB-CAT-BDD-TDR-ITT-11

**DISPOSABLE INTRAVENOUS
CATHETER**

INFUSION / TRANSFUSION
THERAPY

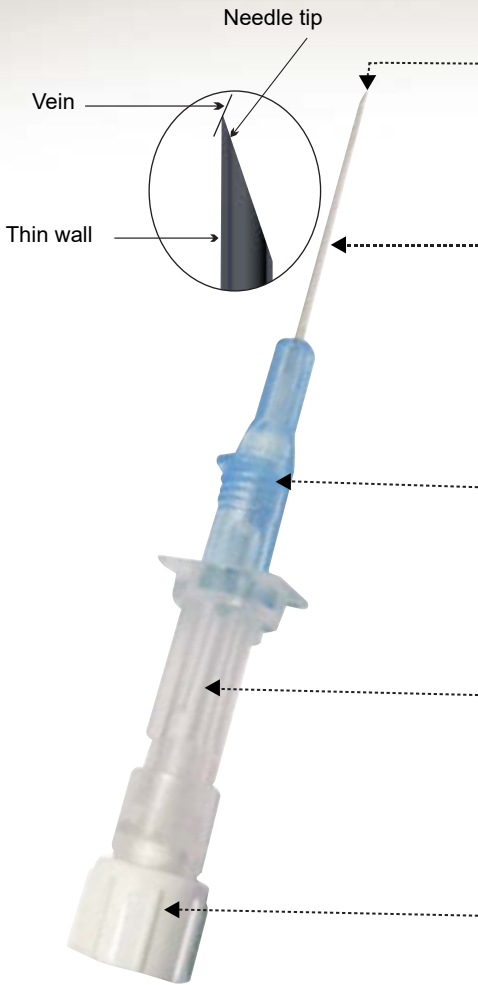

Sharp 3-Facet Back-Cut Needle

- Sharp and thin walled needle with reduced bevel angle for less contact area during insertion resulting **atraumatic cannulation**.

Biocompatible Polyurethane (PUR) Catheter

- Special tapered bio-compatible **Polyurethane (PUR)** catheter reduces the chances of phlebitis and thrombo-phlebitis; kink resistant allow longer indwelling time.
- Smooth surface catheter embedded with **radiopaque lines** for better x-ray visualization.

Catheter Hub

- Ergonomic product design with finger guard facilitates easy insertion **with one-handed insertion**.
- Universal color-coded on catheter hub for easy gauge size identification.

Transparent Flashback Chamber and Hydrophobic Filter Membrane

- Transparent flashback chamber for **easy visualization** of blood.
- Hydrophobic filter to prevent blood leakage during cannulation.
- Ribbed design on chamber secure handling during cannulation.

Removable Closing Cone

- Luer lock plug at the end of catheter to **prevent blood contamination**.

- ✓ Disposable.
- ✓ Sterile.
- ✓ Easy to use.
- ✓ Latex free.
- ✓ PVC free.
- ✓ Individually packed in hard blister packaging.
- ✓ Packaging: 50pcs/box, 10 boxes / carton.

I.V. CATHETER WITHOUT PORT AND WITHOUT WINGS	CATALOGUE NO.	SIZE	COLOR CODE	OD (mm)	ID (mm)	LENGTH (mm)	FLOW RATE (ml/min)
	IC2003-14	14G	Orange	2.0	1.7	45	305
	IC2003-16	16G	Grey	1.7	1.3	45	200
	IC2003-18	18G	Green	1.3	0.9	45/32	90
	IC2003-20	20G	Pink	1.1	0.8	32	56
	IC2003-22	22G	Blue	0.9	0.6	25	36
	IC2003-24	24G	Yellow	0.7	0.5	19	23
	IC2003-26	26G	Purple	0.6	0.45	19	17



Ideal Healthcare Sdn. Bhd. (484861-V)
No. 70 & 71, Jalan Sungai Tukang 2/1, Kawasan Perusahaan Sungai Tukang,
08000 Sungai Petani, Kedah, Malaysia.
Tel : 6 (04) - 422 7775 / 422 8775 Fax : 6 (04) - 422 6775
Email : info@idealcare.com.my

E. L. NO.: MDA-0003-KD314

www.idealcare.com.my

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= ISO 13485 =



idealcareMY

IHSB-CAT-BDD-TDR-ITT-12